

Specifications 220V-240V/50Hz, 5A Power supply 350W Working power 300W Heating power Max.Temp. 200℃ Heating mode Block heating Cooling mode Water cooling 0.5L-2L Feed capacity HNO₃ 1.8 L; HCL 1.6 L; HF 1.5 L; H₂O 2.0 L Production rate (24 hours) Receiving vessel 500mL Distillation rate 1.25mL/min 361 mm(H) × 360mm(W) × 400mm(D) Dimensions

Subject to change without notice. PreeKem reserves rights of final explanation. Version: 1st edition, 2018.12

PreeKem Scientific Instruments Co., Ltd.



P.C.: 201108

Tel: +86-21-54427296 54426316 54426318

Fax: +86-21-54427063

E_mail: info@preekem.com

Web: www.preekem.com





Sub-boilng Acid Purifier

S1 Sub-boiling Acid Purifier


S1 sub-boiling acid purifier adopts classic sub-boiling distillation approach and works to achieve a dynamically delicate balance between automatic control of heat convection and acid vaporization.

The instrument control is designed to enable intelligent operating and handling of any task, making even an unmanned acid purification session possible in the fume hood. The device is constructed with absolutely zero quartz or metal parts, entirely eliminating the occurrence of corrosion and outside contamination from happening.

Performance

3D Surrounding heating mode

3D surrounding heating mode to improve energy utilization in a more efficient and even manner.

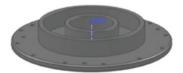

Direct acid injection

Feed acid to be injected from the bottom via specially designed delivery channel to reduce length of passage and eliminate contamination.

Submerged temperature monitoring

Pt100 temp sensor to provide accurate control of acid heating.

Wide-angle spherical cooling chamber

Optimized geometry of cooling chamber to provide highly efficient condensation.



Intuitive control software

Built-in distillation method archive, real-time display of distillation progress.



S1 Sub-boiling Acid Purifier

6

Enhanced safety features



Automatic temperature regulation

Disabling heating when a pre-set temperature threshold is reached and restarting heating when temperature falls back in normal working range



Over heating prevention

System stops automatically if feed acid is below safe working volume.



Intelligent pressure release gauge

Release excess vapor when the internal pressure alarm is triggered.



Strong carrying handle

Providing handing and operating an extra layer of protection.



Direct view of acid level

keeping operators in the know of the distillation process.



Steady base of receiving vessel specially designed to provide extra stability of the receiving vessel while in working mode.



Unparalleled rate of production!



HCL 1.6 L





1 2